

ESSENTIALLY BETTER




Mobilo

The stretcher for flexible use



▲ As a stretcher and day bed, the Mobilo offers patients particularly high safety and comfort. This is ensured by features such as the 35 cm high safety sides on either side of the stretcher and the comfortable cold foam mattress.

External dimensions/height adjustment range:

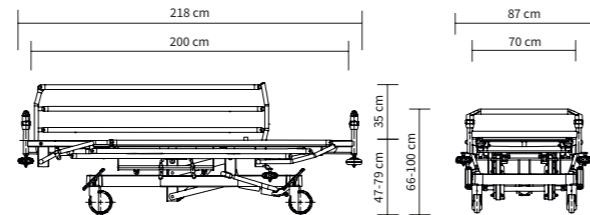

▲ The mattress base height and position can be flexibly adjusted from both sides - from the treatment position right up to a comfortable sitting position.



▲ Safety sides present additional relief for nursing staff as, depending on the situation, the patient no longer requires constant supervision.



▲ The mattress base height can be adjusted quickly and easily using a pedal.



▲ In an emergency, the Mobilo can be tilted swiftly to a Trendelenburg position.

The Mobilo

More than meets the standard

The Mobilo Stretcher serves many important purposes at once. As a **stretcher for transporting** patients to functional areas, it is designed to move patients swiftly and safely from one place to another. With its extensive features, it additionally meets standard DIN EN 60601-2-52 for beds. It can also be **used as a bed in day-care centres**, if required. This provides patients with a comfortable positioning solution. Depending on their diagnosis and state of health, patients are able to spend prolonged periods in the Mobilo without constant supervision.

In addition, the Mobilo facilitates the work of nursing staff. Complex reprocessing that is needed for beds requiring a linen change is not necessary with the stretcher. Its **central braking system** can be operated with brake pedals on all 4 castors. Particularly in stressful situations, nursing staff can brake or move the stretcher from any direction.

Safety and comfort

Support for patients and care staff

The Mobilo mattress base is available in a two-section or four-section version. With its 70 x 200 cm, it is very generously sized for a stretcher. With its **height adjustment** range of 47 to 79 cm, patients are assured of a comfortable entry and exit height, while nursing staff appreciate the ergonomic working height. The swivel-down **safety sides** offer the patient a high level of safety.

In an emergency, when every second counts, the optional **resuscitation support** simply folds out from under the mattress base at the head end. During resuscitation, it provides a stable, vibration-free base.



The tray for resting on the safety sides facilitates patient care without the need for repositioning



▲ The wide, ergonomically designed push bars at the head and foot end make it easy to manoeuvre the stretcher swiftly and easily.

Care in every situation with useful accessories

All care items required for treatment can be securely attached to the Mobilo both quickly and easily. This allows the Mobilo to be used flexibly as a stretcher or day bed.



▲ The optional infusion stands are of double telescopic design and can be fitted to all four corners of the stretcher.

▼ Essential care accessories can be versatilely attached to the supply rails on both sides of the Mobilo.



▼ A storage tray on the chassis offers an easily accessible space for holding an oxygen bottle with its strap and the patient's personal belongings.



Extras for many requirements

With a versatile range of accessories and numerous details, the stretcher can be adjusted to many demands of the daily hospital routine. The easy **manoeuvrability** and excellent **hygiene**

properties of the Mobilo relieve the strain on staff and save time that can be used instead to care for the patient.

▼ The paper roll holder can be fitted either at the head or foot end and ensures particularly hygienic use of the mattress.



▼ Supply accessories can be attached to the sliding rails at the head and foot end of the stretcher.



▼ Thanks to its radiolucent backrest and optional cassette holder, it is no longer necessary to reposition the patient.



Technical data

Equipment and options	
Safety side	
With safety sides on both sides	✓
Basic colour	
Argentum similar to RAL 9006	✓
Mattress bases	
2-section mattress base	✓
4-section mattress base	0
Wall deflection rollers	
6, horizontally and vertically effective	✓
Castors with central locking	
Double castor Ø 150 mm	✓
Integral castor Ø 200 mm	0
Directional castor positions	
Head end, right	✓
Head end, left	0
Foot end, right	0
Foot end, left	0
Conductive castor coating	
1 castor	✓
2 castors	0
4 castors	0
Additional features	
Stretcher tray and X-ray cassette holder under the mattress base	0
Resuscitation support	0

Technical data	
Dimensions and weights	
Overall dimensions	87 x 218.5 cm
Mattress base (mattress dimensions)	70 x 200 cm
Mattress height	10 cm
Safe working load	200 kg
Maximum lifting capacity	250 kg
Ground clearance with mattress base in lowest position	approx. 9 cm
Total weight	approx. 130 kg
Mattress base in 2 sections	74 cm – 117.5 cm
Mattress base in 4 sections	74 cm – 13 cm – 37.7cm – 61.8 cm
Protection height of safety sides	35 cm
Adjustment ranges	
Tilting to Trendelenburg position	12°
Tilting to reverse-Trendelenburg position	12°
Mattress base height	approx. 47 – 79 cm
Backrest angle	70°
Thigh-rest angle (with mattress base in 4 sections)	45°



Custodians of values Shapers of the future

STIEGELMEYERGROUP

Custodians of values

Our company was founded in Herford in 1900 and is still a family-run enterprise today. Respect and a sense of responsibility towards our customers, employees and partners form the basis for our day-to-day work. Based on these values and our many years of experience, we offer our customers all over the world high-quality products, solutions and services to meet the needs of the care and healthcare sectors.

Shapers of the future

In order to cope with rapid changes in the care and healthcare sectors, we develop all our products ourselves. Creating innovations and visions of what care will be like in future are the goals we pursue. We maintain regular contact with our customers in order to support them in their daily tasks. Implementing tomorrow's trends today is a successful tradition of ours.

Our expertise

Hospital

Nursing home

Homecare

Lifestyle beds

Furniture systems

Service



The Stieglmeyer-Group is your international partner for hospitals and care facilities and, with the Burmeier brand, also for home care. We support our customers with practical products, solutions and services, which we constantly develop further.

- Founded in 1900 and still family-run today
- Employs 1,100 members of staff worldwide
- Production of all core products in our own factories
- Product development and manufacture in Germany
- Services, training and process optimisation at the customer's site
- International certification and export to more than 60 countries
- Sales companies in Belgium, Finland, France, the Netherlands, Poland and South Africa

Our brands



STIEGELMEYERGROUP

Hauptsitz / Headquarters

Deutschland / Germany
Stieglmeyer GmbH & Co. KG
Ackerstraße 42
32051 Herford
Phone +49 (0) 5221 185 - 0
Fax +49 (0) 5221 185 - 252
info@stieglmeyer.com
www.stieglmeyer.com

France
Stieglmeyer SAS
8 Allée des Chevreuils
69380 Lissieu
Phone +33 (0) 4 7238 - 5336
Fax +33 (0) 4 7238 - 5481
contact@stieglmeyer.fr
www.stieglmeyer.fr

South Africa
Stieglmeyer Africa (PTY) Ltd
The Estuaries · Century City
Cape Town, 7441
Phone +27 (0) 21 555 - 1415
Fax +27 (0) 21 555 - 0639
info.southafrica@stieglmeyer.com
www.stieglmeyer.co.za

Nederland / Netherlands
Stieglmeyer BV
Peppelenbos 7
6662 WB Elst
Phone +31 (0) 345 581 - 881
info@stieglmeyer.nl
www.stieglmeyer.nl

Suomi / Finland
Stieglmeyer Oy
Ahjokatu 4 A
15800 Lahti
Phone +358 (0) 3 7536 - 320
Fax +358 (0) 3 7536 - 330
stieglmeyer@stieglmeyer.fi
www.stieglmeyer.fi

Polska / Poland
Stieglmeyer Sp. z o.o.
86-212 Stolno
Grubno 63
Phone +48 56 6771 - 400
Fax +48 56 6771 - 410
info.stolno@stieglmeyer.pl
www.stieglmeyer.pl



Finden Sie Ihren Ansprechpartner auf unserer Website:
Find your local partner on our website:
www.stieglmeyer.com/partner

stieglmeyer.com