



Classiko

Robust, yet homely



Classiko

The comfortable and secure solution

Comfort and safety for the resident, relief for caregivers, and an individual design – the Classiko unites many strengths to give good value for money. This robust care bed has the advantage of durability and a range of adjustment options that make everyday life more agreeable.

The four-section mattress base adapts to suit many different needs. For example, by raising the backrest and thigh rest and tilting the bed into the optional reverse-Trendelenburg position, the resident is moved into a comfortable sitting position for eating or watching television. With a height adjustment range from 41 to

78 cm, the Classiko makes getting in and out of bed easy while ensuring a back-friendly working position for care staff. The electrical model of the bed can be operated effortlessly with a selectively lockable handset.

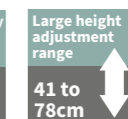
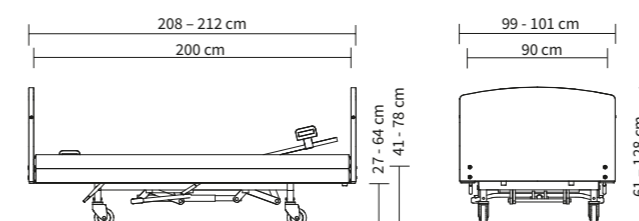
An electrically adjustable bed is not necessary or desirable in all care situations, however. For this reason, Stiegmeyer also offers a mechanical version of the Classiko. By dispensing with cables and manually removable elements, this model meets the safety requirements in psychiatric care particularly well. With a high safe working load of 225 kg, it is also suitable for heavy residents.

Care staff will find the mechanical Classiko easy to operate. Height adjustments using a pedal with the help of a hydraulic pump can be made just as easily as the power-assisted adjustments of the backrest and thigh rest. The operator makes savings in energy costs and expenses for leakage current measurements. The Classiko offers visual appeal with a large selection of different headboards and footboards and wood decors to suit every taste.

ESSENTIALLY BETTER
High safe working load

- electrically adjustable: 175 kg
- mechanically adjustable: 225 kg

External dimensions/height adjustment range:





The Classiko bed is fitted with a full-length safety side as standard.

Optimum protection for the resident

Safety sides and practical equipment options

The safety sides on the Classiko can be adjusted to meet the individual protection needs of the resident. The standard version of the bed is equipped with full-length safety sides on both sides of the mattress base for optimum safety.

If the resident's mobility needs to be supported, then the combinable split safety sides are a good choice. When raised at the head end only, it serves as a stable support for getting in and out of bed. The split safety side can also be combined with a full-length safety side on the other side of the bed.

When safety sides are not required at all – the Classiko can be fitted with a fixed side panel, which gives the bed a particularly homelike appearance.

With its features, options and accessories, the Classiko becomes a comfortable bed adapted to your requirements. The four-section steel wire mesh mattress base is robust and flexible. In the electrically model, the auto contour setting and the auto-regression on the backrest minimise shear.

For taller residents, we recommend bed extensions to 210 or 220 cm. Matching mattress extension pieces and extended safety side bars are also available.

The Classiko is easy to move and manoeuvre. Thanks to its central-locking castors, care staff can lock or unlock all four castors at once by pressing just one brake pedal. The electrically model is operated using a handset, which shows all of the options as easily understandable pictograms. Individual functions can be blocked for the resident's safety, if required.



▲ Without safety sides the Classiko is ideal for independent residents.



▲ The 4-section wire mesh mattress base is both sturdy and flexible.

▼ The Classiko care bed is optionally available with combinable split safety sides.



▼ The bed extension provides extra comfort for tall residents.



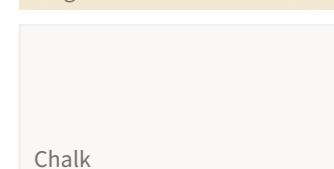
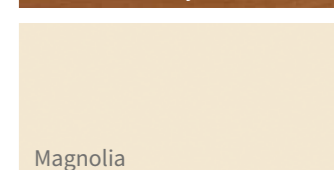
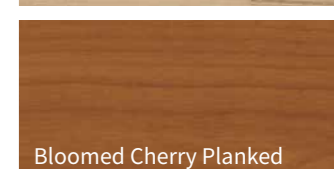
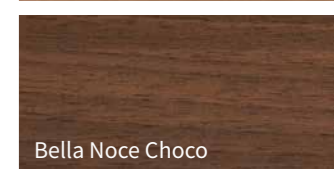
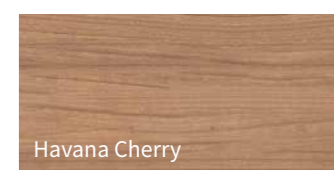
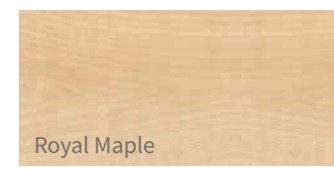
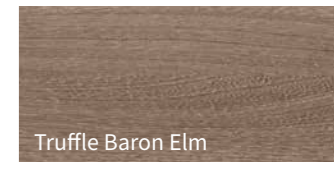
▲ Bed adjustments on the mechanical model use a pedal and are hydraulically assisted.





▲ Selectively lockable handset for easier operation.

Decors

The Classiko can be adapted to suit any room with a variety of attractive decors to choose from. The selection ranges from the light Lindberg Oak to the striking dark Bella Noce Choco.



Equipment and options

Dimensions and weights	Mechanical	Electric
External dimensions (depending on the wooden bed frame - W x L)	approx. 101 x 208 cm	approx. 101 x 208 cm
Mattress base (mattress dimensions - W x L)	90 x 200 cm	90 x 200 cm
Backrest length compensation	approx. 10 cm	approx. 10 cm
Safe working load	225 kg	175 kg
Ground clearance	approx. 15 cm	approx. 15 cm
Total weight	approx. 130 kg	approx. 130 kg
Standard features		
4-section mattress base	✓	✓
Wire mesh mattress base	✓	✓
Sheet steel castor Ø 125 mm	✓	✓
Directional castor, head end, right	✓	✓
4 central locking castors	✓	✓
Height adjustment range 41 to 78 cm	✓	✓
Backrest adjustment to 70°	✓	✓
Extended leg raise	✓	✓
Auto contour setting	-	✓
Handset with selective locking function	-	✓
Safety side systems		
No safety side (panels only)	✓	✓
Full-length safety side 	0	0
Combinable safety side 	0	0
Accessory equipment		
Trendelenburg and reverse-Trendelenburg position	0	0
Bed extension to 210 or 220 cm	0	0
Universal accessory holder	0	0
Secured pump lever	0	-
Attachment points for posey belts on mattress frame and/or mattress base	0	0
Wall deflection roller, all sides - vertical or universal	0	0
Wall deflection roller, head end - vertical or universal	0	0
Wall deflection roller, universal, at head end and vertical, at foot end	0	0

✓ = Standard feature
 0 = Available as an option
 - = Not available

Safety and comfort for psychiatric care

The mechanical model of the Classiko is perfectly suited to the high safety requirements in psychiatric care. Since it has no cables and its components are protected against removal, it prevents injuries to residents and damage to the bed. At the same time, the comfort and homelike appearance of the Classiko enhance the feeling of well-being.

For more information on use in psychiatric care, go to www.stieglmeyer.com

STIEGELMEYERGROUP



Solutions for people who care

The Stieglmeyer Group is Germany's leading manufacturer of medical beds for hospitals, nursing homes and homecare. With modern bedside cabinets, furniture and digital assistance systems, it is also an innovative provider of comprehensive solutions. The family-run company successfully develops offers for international markets in more than 100 countries.

Our strengths:

- Product development at our Herford headquarters in Germany – production in four modern factories in Germany and Poland
- Factory in China for the Chinese market
- Subsidiaries and sales partners in many countries
- Top-class service with reliable spare parts supply
- Training and process optimisation for our customers
- Sustainable action with strong social commitment

1,000 EMPLOYEES WORLDWIDE

11 LOCATIONS in Germany, Poland, the Netherlands, Belgium, France, Finland, South Africa and China

Tailor-made OFFERS for all MARKETS

MORE THAN 100,000 beds produced every year

UP TO 15-YEAR spare part availability

DISTRIBUTION in more than 100 countries

INTERNATIONAL CERTIFICATION

GERMAN ENGINEERING

Our expertise

- Hospital
- Nursing home
- Homecare
- Lifestyle beds
- Furniture systems
- Service

Our values

- Experienced
- Human
- Reliable
- Solution-focused
- Open-minded
- Confident

Our brands





STIEGELMEYERGROUP

Solutions for people who care



Hauptsitz / Headquarters

Deutschland / Germany
Stieglmeyer GmbH & Co. KG
Ackerstraße 42
32051 Herford
Phone +49 (0) 5221 185-0
Fax +49 (0) 5221 185-100
info@stieglmeyer.com
www.stieglmeyer.com



Hauptsitz / Headquarters

Deutschland / Germany
Burmeier GmbH & Co. KG
Industriestraße 53
32120 Hiddenhausen
Phone +49 (0) 5223 9769 - 0
Fax +49 (0) 5223 9769 - 090
info@burmeier.com
www.burmeier.com

France
Stieglmeyer SAS
16 Avenue de l'Europe
Bâtiment SXB1
67300 Schiltigheim
Phone +33 (0) 4 7238 - 5336
contact@stieglmeyer.fr
www.stieglmeyer.fr

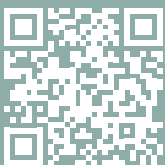
Nederland / Netherlands
Stieglmeyer BV
Peppelenbos 7
6662 WB Elst
Phone +31 (0) 345 581 - 881
info@stieglmeyer.nl
www.stieglmeyer.nl

Polska / Poland
Stieglmeyer Sp. z o.o.
86-212 Stolno
Grubno 63
Phone +48 56 6771 - 400
Fax +48 56 6771 - 410
info.stolno@stieglmeyer.pl
www.stieglmeyer.pl

South Africa
Stieglmeyer Africa (PTY) Ltd.
The Estuaries - Century City
Cape Town, 7441
Phone +27 (0) 21 555 - 1415
Fax +27 (0) 21 555 - 0639
info.southafrica@stieglmeyer.com
www.stieglmeyer.co.za

Suomi / Finland
Stieglmeyer Oy
Ahjokatu 4 A
15800 Lahti
Phone +358 (0) 3 7536 - 320
Fax +358 (0) 3 7536 - 330
stieglmeyer@stieglmeyer.fi
www.stieglmeyer.fi

China
Stieglmeyer China Co. Ltd.
Poly Center, No. 8 South Wenhua Road
Dongping New Town, Shunde District
528300 Foshan
Phone +86 757 29878710 - 0
info@stieglmeyer.cn
www.stieglmeyer.cn



stieglmeyer-group.com