



Carino

Feel at home wherever you are



▲ The Carino in a combination of a Nobla reha headboard and a Nobla mala footboard, with a side tray attached to the bed frame.

Carino

Fresh ideas for comfort and stability

Stable, durable, comfortable and extremely versatile – the Carino is a bed with many possible applications. It greatly enhances the quality of life in rehabilitation facilities as well as in psychosomatic treatment centres, residential homes and initial accommodation for migrants. The Carino is particularly appreciated for its very good price-performance ratio and individual design possibilities.

With its robust metal mattress base, the bed easily withstands numerous cycles of use. Mattress base widths of 90 and 100 cm are available and a bed extension can be supplied for taller users. The Carino's manually adjustable backrest ensures a comfortable position when reading or watching TV. The optional reading lamp provides pleasant lighting – fixing sleeves are already in place and can also be used for fixing a patient lifting pole. Another practical accessory is the space-saving side tray, providing convenient space for a book or mobile phone.

The Carino has a sleek and elegant design and can be adapted to your own personal taste with a wide range of available headboard and footboard forms and wood decors to choose from. The textile Softcovers for the headboard, footboard and sides are an especially cosy design feature of the bed. Fabrics in attractive colours and patterns can be used to transform the Carino into the centrepiece of the room. With an optional fixed board along the side of the bed, the Carino can also function as a comfortable sofa during the day.



Optional bed drawers offer storage space under the bed. ▶

Carino	
Dimensions and weights	
External dimensions (depending on the wooden bed frame)	approx. 105 x 210 cm or 115 x 210 cm
Mattress base (mattress dimensions)	approx. 90 x 200 cm or 100 x 200 cm
Mattress base height	approx. 40 cm
Mattress base division (backrest, mattress base)	approx. 80/120 cm
Safe working load	225 kg
Ground clearance	25 cm
Total weight *	approx. 100 kg
Adjustment ranges	
2-section mattress base, manually adjustable backrest	✓
Backrest adjustment to 30°	✓
Miscellaneous	
Support elements with easy-to-clean metal slats	✓
2 adapter sleeves at the head end for patient lifting pole/reading lamps	✓
Wooden surrounds, choice of Lito / Lito reha, Alero, Nobla / Nobla reha / Nobla mala, Pino, Podego, Pilastro, (dimensions vary depending on the model of headboard and footboard)	✓
Choice of wood decors from Stieglmeyer Wodego collection	✓
Coated metal finish in Iron Grey RAL 7011	✓
Feet, chrome-plated	✓
Accessories	
Side tray	0
Softcovers	0
Bed drawer	0
Reading lamp, LED light bulb	0
Additional features, accessory equipment	
Staining of wooden surround	0
Mattress base width 100 cm	0
Bed extension 20 cm **	0
Model with fixed board along one side of the bed	0
Selected headboards/footboards in multiplex material	0
Side panels, attached with through bolts	0
Patient lifting pole with grab handle	0
Cable holder for patient lifting pole	0
Infusion stand, straight or angled	0

Particularly sturdy mattress base and 225 kg safe working load

✓ = Standard feature
 0 = Available as an option
 - = Not available

* Depending on the materials and accessories used, the overall weight may vary from the figure shown.

** Not with 100 cm mattress base



▼ The metal slats of the mattress base are easy to clean.

◀ The fixed board along the side of the bed transforms the Carino into a sofa during the day.



STIEGELMEYERGROUP

Solutions for people who care



Hauptsitz / Headquarters

Deutschland / Germany
Stieglmeyer GmbH & Co. KG
Ackerstraße 42
32051 Herford
Phone +49 (0) 5221 185-0
Fax +49 (0) 5221 185-100
info@stieglmeyer.com
www.stieglmeyer.com



Hauptsitz / Headquarters

Deutschland / Germany
Burmeier GmbH & Co. KG
Industriestraße 53
32120 Hiddenhausen
Phone +49 (0) 5223 9769 - 0
Fax +49 (0) 5223 9769 - 090
info@burmeier.com
www.burmeier.com

France

Stieglmeyer SAS
16 Avenue de l'Europe
Bâtiment SXB1
67300 Schiltigheim
Phone +33 (0) 4 7238 - 5336
contact@stieglmeyer.fr
www.stieglmeyer.fr

Nederland / Netherlands

Stieglmeyer BV
Peppelenbos 7
6662 WB Elst
Phone +31 (0) 345 581 - 881
info@stieglmeyer.nl
www.stieglmeyer.nl

Polska / Poland

Stieglmeyer Sp. z o.o.
86-212 Stolno
Grubno 63
Phone +48 56 6771 - 400
Fax +48 56 6771 - 410
info.stolno@stieglmeyer.pl
www.stieglmeyer.pl

South Africa


Stieglmeyer Africa (PTY) Ltd.
The Estuaries - Century City
Cape Town, 7441
Phone +27 (0) 21 555 - 1415
Fax +27 (0) 21 555 - 0639
info.southafrica@stieglmeyer.com
www.stieglmeyer.co.za

Suomi / Finland

Stieglmeyer Oy
Ahjokatu 4 A
15800 Lahti
Phone +358 (0) 3 7536 - 320
Fax +358 (0) 3 7536 - 330
stieglmeyer@stieglmeyer.fi
www.stieglmeyer.fi

China

Stieglmeyer China Co. Ltd.
Poly Center, No. 8 South Wenhua Road
Dongping New Town, Shunde District
528300 Foshan
Phone +86 757 29878710 - 0
info@stieglmeyer.cn
www.stieglmeyer.cn

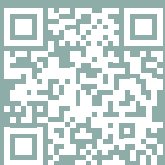

stieglmeyer-group.com