



The care bed with innovative ideas


Inovia

ADVANTAGES

The Inovia care bed has an attractive, homelike appearance thanks, in particular, to its opulent headboards and footboards. Its lowest position, at just 22 cm, protects resident from falls, without having to resort to liberty-depriving measures. The optional telescopic safety sides also promote the resident's mobility.

- Safe working load: 225 kg
- Maximum weight of resident: 195 kg
- Height adjustment range: 22 to 77 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Reverse-Trendelenburg position
- Sliding backrest in accordance with the recommendations of the German Nurses Association (DBfK)
- Optional telescopic split safety side
- Attractive wood decors to choose from

External dimensions/height adjustment range:

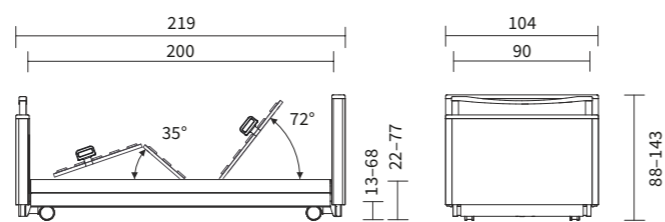

Special features:



Further information:



The advantages at a glance

DESIGN DETAILS



Solid wood handle bar.



Castors can be locked in pairs whatever the mattress base position.



Handset with selective locking.



Sliding backrest.



Mattress base height adjustable from 22 to 77 cm, adjustment range = 55 cm.



Natural Beech (standard).



Havana Cherry.

ACCESSORIES



Telescopic safety side (TSG).



Foam leather cover for DSG.



Foam leather cover for TSG.



Bed extension.



Tray.



Sola reading lamp.



The storage tray for glasses, mobiles or care utensils can be attached to the full-length safety side.



The comfort mattress base comprises 50 free-floating spring elements.



Flexible transmitter holder for the handset.

EQUIPMENT AND OPTIONS

Equipment/Option	Inovia	Material no.
Design features		
Inovia 90 x 200 cm with full-length safety sides (DSG)		262626LFB_001 SHC_001
Inovia 100 x 200 cm with full-length safety sides (DSG)		262626LFB_002 SHC_001
Technical data		
Net weight	153/159 kg	
Overall dimensions	105 x 220 cm 115 x 220 cm	
Height adjustment range	22 – 77 cm	
Safe working load	225 kg	
Maximum resident weight	195 kg	
24-volt system	✓	
Standard features		
Frame colour	Argentum	
Decor	Natural Beech	
Mattress base size	90 x 200 cm 100 x 200 cm	
Metal slatted mattress base (mattress base model)	✓	
Type of handset (selectively lockable/wireless)	Selective	
Mechanical CPR release	✓	
Castor diameter	100 mm	
Castor locking	Per axle	
Ratchet mechanism, electrical adjustment of lower leg rest	✓	
Mechanical CPR release	✓	
Special equipment and accessories		
Telescopic safety sides (TSG) on both sides	0	SHC_002
Telescopic safety side (TSG) on one side/DSG on one side	0	SHC_008
Patient lifting pole	0	AUI_001
Natural Beech	✓	DHU_027
Havana Cherry	0	DHU_024
Handle bar at head end	0	KTA_003
Bed extension for metal slatted mattress base 90 and 100 cm with TSG	0	262626LFL_006
Bed extension for metal slatted mattress base 90 and 100 cm with DSG	0	262626LFL_005
Mattress piece for bed extension 90 cm wide, 12 cm thick	0	166715
Mattress piece for bed extension 100 cm wide, 12 cm thick	0	219857
Comfort mattress base	0	LIF_013
Wall deflection roller with holder (right or left)	0	Left 215921 Right 221306
Foam leather cover for full-length safety sides (DSG)	0	225921
Foam leather cover for telescopic safety sides (TSG)	0	716556
Reading lamp for easy insertion	0	264877
Reading lamp for fixing to the patient lifting pole	0	1 x 191879 1 x 174231
Tray for full-length safety sides (DSG), 90 cm wide	0	174296
Tray for full-length safety sides (DSG), 100 cm wide	0	229222
Extension cushion (for full-length safety side only)	0	258395
Transmitter holder for handset	0	199137
Handset extension cable	0	238047
Storage tray	White Yellow Red	266970 266971 266619

✓ = Standard feature 0 = Available as an option