



The innovative care beds


Allura II

ADVANTAGES


The flexible models of the Allura II range all have a high load capacity with a choice of mattress base widths and a large height adjustment range. This means that they offer even heavy residents a high level of comfort and safety, while relieving the burden placed on carers.

- Safe working load: 250 kg
- Maximum weight of resident: 210 kg
- Height adjustment range: 30 to 80 cm
- 24-volt drive system
- No leakage current measurements during expected service life
- Reverse-Trendelenburg position
- Selectively lockable handset
- Mattress base widths of 100 and 120 cm
- Tool-free safety side installation with Easy Click system

External dimensions/height adjustment range:

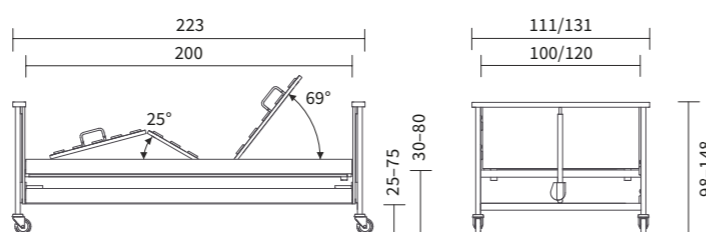

Special features:

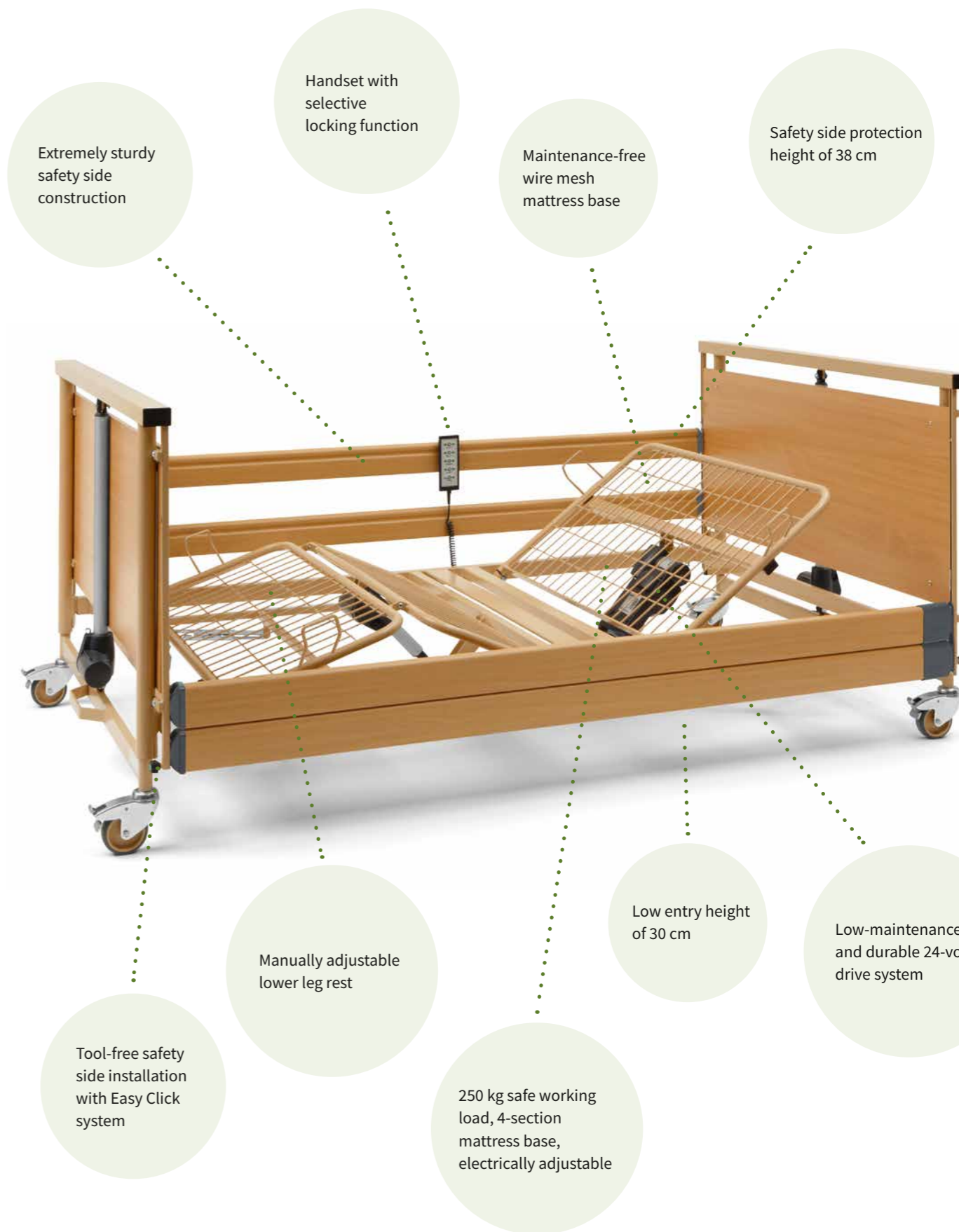

Further information:



The advantages at a glance

DESIGN DETAILS



Reverse-Trendelenburg position.



Adjustment range from 30 to 80 cm.



Electronic (SwitchMode) power supply unit integrated immediately downstream of the mains plug.



Handset with selective locking function.



Compact delivery and storage on storage aid/transportation aid (similar to illustration).



Tool-free safety side installation with Easy Click system.

ACCESSORIES



Bed extension (also for retrofitting).



Foam leather cover.



Reading lamps.



Wooden panels to cover the chassis (also for retrofitting).



Extension cushion to place over the full-length safety side.



Wall deflection roller including holder.



Standing aid.

EQUIPMENT AND OPTIONS

| Equipment/Option | Allura II | Material no. |
|--|-------------------------------|--------------|
| Design features | | |
| Allura II 100x200 cm | | 289450 |
| Allura II 120x200 cm | | 289451 |
| Technical data/standard features | | |
| Net weight | 120 kg | |
| Overall dimensions | 111 x 224 cm | |
| Height adjustment range | 30 – 80 cm | |
| Safe working load | 250 kg | |
| Maximum resident weight | 210 kg | |
| 24-volt drive system | ✓ | |
| Standard features | | |
| Frame colour | Topaz | |
| Decor | Natural Beech | |
| Mattress base size | 100 x 200 cm | |
| Safety side | Full-length safety side (DSG) | |
| Type of handset (selectively lockable/wireless) | Selective | |
| Mechanical CPR release | ✓ | |
| Castor diameter | 100 mm | |
| Castor locking | Single | |
| Ratchet mechanism, electrical adjustment of lower leg rest | ✓ | |
| Accessories | | |
| Standing aid | O | 244029 |
| Bed extension for metal slatted mattress base 100 cm | O | 291745 |
| Bed extension for metal slatted mattress base 120 cm | O | 291748 |
| Mattress piece for bed extension 120 cm wide, 16 cm thick | O | 241942 |
| Wall deflection roller, 1 piece, including mount | O | 206751 |
| Foam leather cover for full-length safety sides (DSG) | O | 225921 |
| Reading lamp for easy insertion | O | 264877 |
| Reading lamp for fixing to the patient lifting pole | O | 289177 |
| Extension cushion (for full-length safety side only) | O | 258395 |