

Out-of-Bed flex

The accommodating solution
for bed exit notifications




Out-of-Bed flex – the sensor bar for flexible use

The Out-of-Bed flex system enables hospitals to cater appropriately to the needs of fall-prone patients. The retrofittable OOB flex assistance system notifies nurses whenever a patient gets up and does not return to bed within a specified period of time. It offers protection and quality of life combined with cost-efficient flexibility in hospital care.

The system can be installed, only when required, on different Stiegmeyer hospital beds. The OOB flex system enables hospitals to meet a tailored investment need, and provide a high standard of care for their patients while also easing the workload of nurses – particularly in situations with fewer people on duty.

Easy to use and install

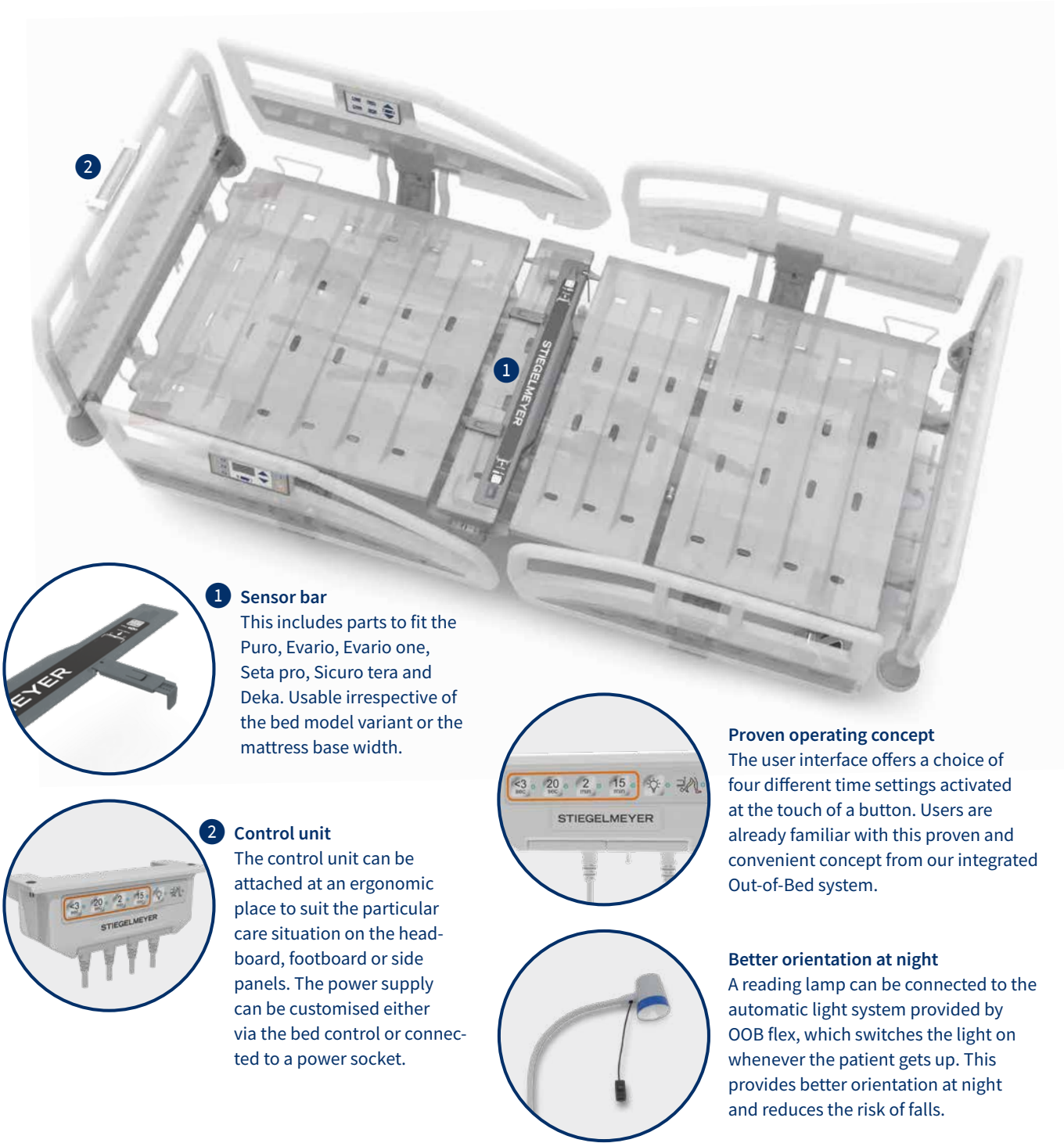

1 Sensor bar
This includes parts to fit the Puro, Evario, Evario one, Seta pro, Sicuro tera and Deko. Usable irrespective of the bed model variant or the mattress base width.

2 Control unit
The control unit can be attached at an ergonomic place to suit the particular care situation on the headboard, footboard or side panels. The power supply can be customised either via the bed control or connected to a power socket.

Proven operating concept
The user interface offers a choice of four different time settings activated at the touch of a button. Users are already familiar with this proven and convenient concept from our integrated Out-of-Bed system.

Better orientation at night
A reading lamp can be connected to the automatic light system provided by OOB flex, which switches the light on whenever the patient gets up. This provides better orientation at night and reduces the risk of falls.



The Out-of-Bed flex system makes work easier for nurses in many different ways. With its control unit placed in reach at an ergonomic working height, there is no need for nurses to look for it or reach below the bed to find it. It also provides a practical aid to decision-making thanks to the system's preset times. These time delays range between under 3 seconds and 15 minutes until an alert is triggered.

During day-to-day nursing, this alleviates the pressure of responsibility while still providing the patient with full protection.

Activating the OOB flex system is possible even when the bed is at its lowest position during patients' rest periods. This makes work appreciably easier for nurses, especially at night.

Your advantages at a glance

- Flexible use – no design or technical requirements to be fulfilled by the Stiegmeyer hospital bed
- Tool-free assembly
- The range of preset times is based on the proven design of the Stiegmeyer OOB smart handset
- Power supply through the use of an external power supply unit
- OOB flex can be enhanced with the addition of a reading lamp
- Flexible option for connecting to standard nurse call systems



Technical data and compatible hospital bed models

OOB flex components and features

Components	Control unit, sensor bar, power supply unit, reading lamp with direct connection cable (optional) and call system connection cable (optional)
Operating voltage	18-40V DC
Protection category	IPX4

Compatible hospital beds

Puro
Evario
Evario one
Seta pro
Sicuro tera
Deka



Hauptsitz / Headquarters

Deutschland / Germany
Stiegmeyer GmbH & Co. KG
Ackerstraße 42
32051 Herford
Phone +49 (0) 5221 185 - 0
Fax +49 (0) 5221 185 - 100
info@stiegmeyer.com
www.stiegmeyer.com



Hauptsitz / Headquarters

Deutschland / Germany
Burmeier GmbH & Co. KG
Industriestraße 53
32120 Hiddenhausen
Phone +49 (0) 5223 9769 - 0
Fax +49 (0) 5223 9769 - 090
info@burmeier.com
www.burmeier.com

France
Stiegmeyer SAS
16 Avenue de l'Europe
Bâtiment SXB1
67300 Schiltigheim
Phone +33 (0) 4 7238 - 5336
contact@stiegmeyer.fr
www.stiegmeyer.fr

Nederland / Netherlands
Stiegmeyer BV
Peppelenbos 7
6662 WB Elst
Phone +31 (0) 345 581 - 881
info@stiegmeyer.nl
www.stiegmeyer.nl

Polska / Poland
Stiegmeyer Sp. z o.o.
86-212 Stolno
Grubno 63
Phone +48 56 6771 - 400
Fax +48 56 6771 - 410
info.stolno@stiegmeyer.pl
www.stiegmeyer.pl

China
Stiegmeyer China Co. Ltd.
Poly Center, No. 8 South Wenhua Road
Dongping New Town, Shunde District
528300 Foshan
Phone +86 757 29878710 - 0
info@stiegmeyer.cn
www.stiegmeyer.cn

South Africa
Stiegmeyer Africa (PTY) Ltd.
5 Shiraz Street, Oude Westhof
Bellville, 7530
Phone +27 (0) 21 555 - 1415
Fax +27 (0) 21 555 - 0639
info.southafrica@stiegmeyer.com
www.stiegmeyer.co.za

Suomi / Finland
Stiegmeyer Oy
Paanakatu 5
15150 Lahti
Phone +358 (0) 3 7536 - 320
stiegmeyer@stiegmeyer.fi
www.stiegmeyer.fi

België / Belgique / Belgium
Distrac Group
Bleyveld 14
3320 Hoegaarden
Phone +32 (0) 16 76 82 82
info@distracgroup.com
www.distracgroup.com



stiegmeyer-group.com